



Anchor for CDS® Hinges

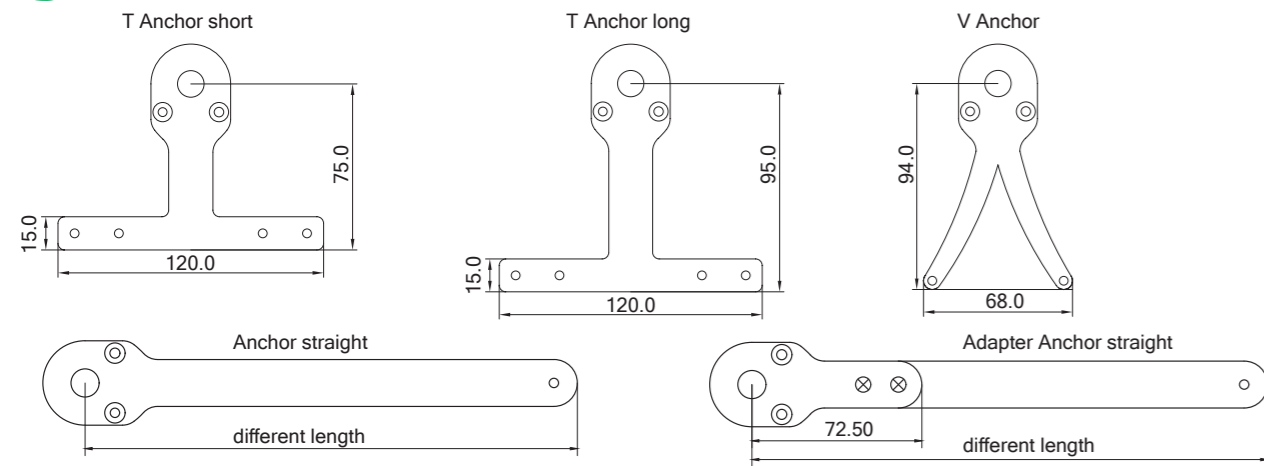


Anchor for CDS® Hinges	Art.No.
CDS® V Anchor stainless steel	813-400-S
CDS® T Anchor short stainless steel	813-405-S
CDS® T Anchor long stainless steel	813-410-S
CDS® V Anchor aluminium	813-420-S
CDS® T Anchor short aluminium	813-425-S
CDS® T Anchor long aluminium	813-430-S
CDS® Anchor straight 16mm stainless steel	813-440-S
CDS® Anchor straight 20mm stainless steel	813-445-S
CDS® Anchor straight 24mm stainless steel	813-450-S
CDS® Anchor straight 20mm aluminium	813-475-S
CDS® Anchor straight 24mm aluminium	813-480-S
CDS® Adapter Anchor straight 16mm stainless steel	813-500-S
CDS® Adapter Anchor straight 20mm stainless steel	813-505-S
CDS® Adapter Anchor straight 24mm stainless steel	813-510-S
CDS® Adapter Anchor straight 20mm aluminium	813-535-S
CDS® Adapter Anchor straight 24mm aluminium	813-540-S

Anchor for CDS® Nano Hinges

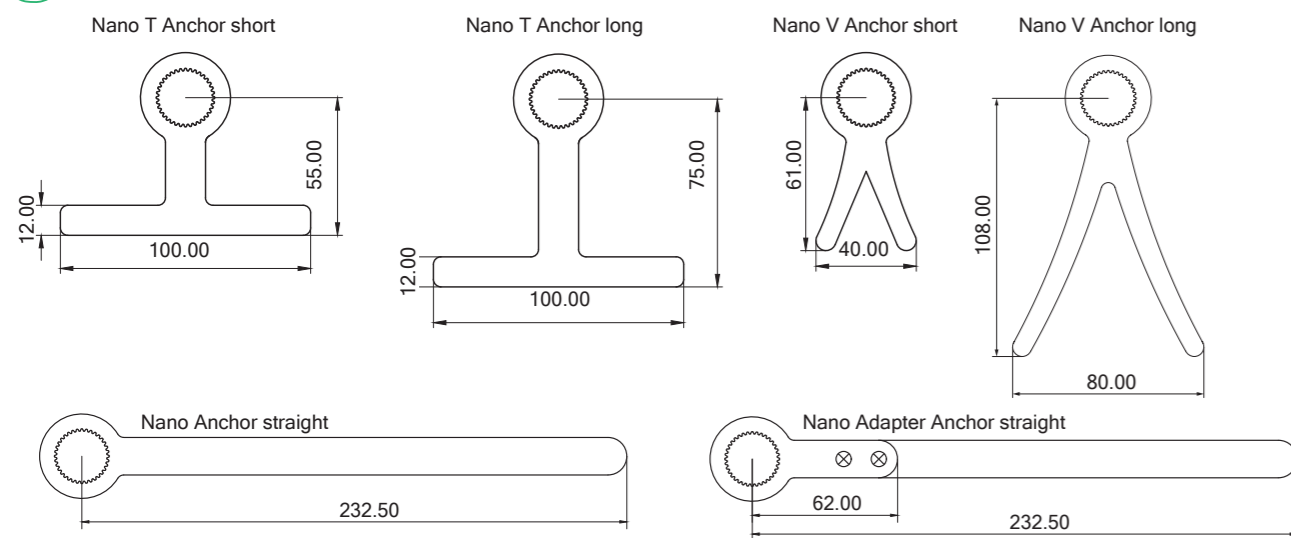

Anchor for CDS® Nano Hinges	Art.No.
CDS® Nano V Anchor short stainless steel	813-100-S
CDS® Nano V Anchor long stainless steel	813-105-S
CDS® Nano T Anchor short stainless steel	813-110-S
CDS® Nano T Anchor long stainless steel	813-115-S
CDS® Nano V Anchor short aluminium	813-130-S
CDS® Nano T Anchor short aluminium	813-135-S
CDS® Nano T Anchor long aluminium	813-140-S
CDS® Nano Anchor straight 16mm stainless steel	813-160-S
CDS® Nano Anchor straight 16mm aluminium	813-170-S
CDS® Nano Adapter Anchor straight 16mm stainless steel	813-190-S
CDS® Nano Adapter Anchor straight 16mm aluminium	813-195-S

CDS® Nano Anchor straight long

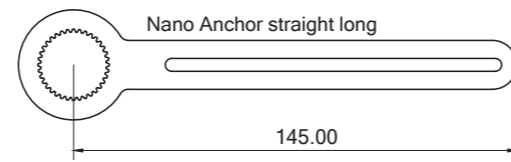

CDS® Nano Anchor straight long	Art.No.
CDS® Nano Anchor straight 16mm long stainless steel	813-165-S
CDS® Nano Anchor straight 16mm long aluminium	813-175-S

Palm Plate Nano



Palm Plate	Art.No. left	Art.No. right
Palm Plate L plus width 10 cm	827LN-PL	827LN-PR
Palm Plate L width 10 cm	827LN-L	827LN-R
Palm Plate M width 8 cm	827MN-L	827MN-R
Palm Plate S width 6 cm	827SN-L	827SN-R

Coniungi



Disconnect Coniungi:
Pull the elements apart while pushing the bolt.



Reconnect Coniungi:
Push the elements into one another.

Coniungi	Art.-No.
Coniungi 16mm	814-100-S
Coniungi 24mm	814-120-S

Stabiliser

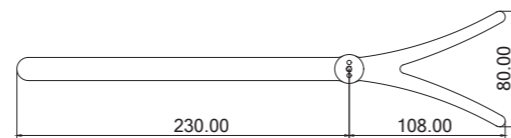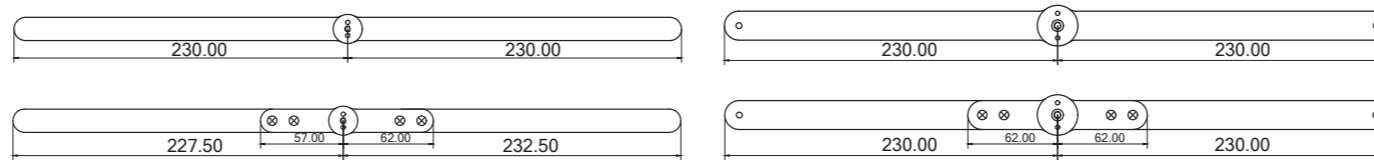

Stabilisers can be configured to match the hinge selected. You can select stabiliser or adapter-stabiliser, aluminium or stainless steel as well as different widths.

Stabiliser	Art.-No.
CDS® Nano Stabiliser stainless steel	830-195-S
CDS® Nano Stabiliser V Anchor long stainless steel	830-196-S
CDS® Nano Stabiliser aluminium	830-190-S
CDS® Nano Adapter Stabiliser stainless steel	830-295-S
CDS® Nano Adapter Stabiliser aluminium	830-290-S
CDS® Stabiliser 16 mm stainless steel	830-180-S
CDS® Stabiliser 20 mm stainless steel	830-160-S
CDS® Stabiliser 24 mm stainless steel	830-110-S
CDS® Stabiliser 20 mm aluminium	830-150-S
CDS® Stabiliser 24 mm aluminium	830-100-S
CDS® Adapter Stabiliser 16 mm stainless steel	830-280-S
CDS® Adapter Stabiliser 20 mm stainless steel	830-260-S
CDS® Adapter Stabiliser 24 mm stainless steel	830-210-S
CDS® Adapter Stabiliser 20 mm aluminium	830-250-S
CDS® Adapter Stabiliser 24 mm aluminium	830-200-S

The CDS® Hinge 360°

Features


- Individual adjustment of the spring tension up to 5,2 Nm
- Spring tension can be switched on and off without tools and without varying the set spring tension
- Range of motion can be limited
- Immobilisation possible

Hinge 360°

- The redression range can be adjusted in 15° increments
- 90° redression range can be positioned freely
- For therapy in extension or flexion
- For treatment of severe extension or flexion deficits

CDS® Hinges 360°	Art.No. left	Art.No. right
CDS® Hinge 360° 16 mm stainless steel	812-200-L	812-200-R
CDS® Hinge 360° 20 mm stainless steel	812-205-L	812-205-R
CDS® Hinge 360° 24 mm stainless steel	812-210-L	812-210-R
CDS® Hinge 360° 20 mm aluminium	812-235-L	812-235-R
CDS® Hinge 360° 24 mm aluminium	812-240-L	812-240-R
CDS® Adapter Hinge 360° 16 mm stainless steel	812-250-L	812-250-R
CDS® Adapter Hinge 360° 20 mm stainless steel	812-255-L	812-255-R
CDS® Adapter Hinge 360° 24 mm stainless steel	812-260-L	812-260-R
CDS® Adapter Hinge 360° 20 mm aluminium	812-285-L	812-285-R
CDS® Adapter Hinge 360° 24 mm aluminium	812-290-L	812-290-R

The CDS® Hinge 360° plus

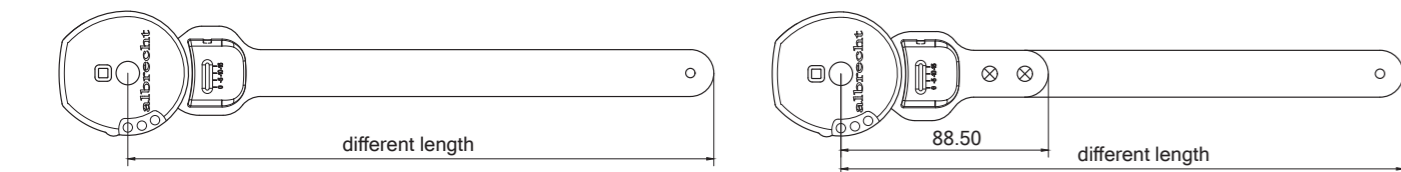

CDS® Hinges 360° plus	Art.No. left	Art.No. right
CDS® Hinge 360° plus 16 mm stainless steel	812-300-L	812-300-R
CDS® Hinge 360° plus 20 mm stainless steel	812-305-L	812-305-R
CDS® Hinge 360° plus 24 mm stainless steel	812-310-L	812-310-R
CDS® Hinge 360° plus 20 mm aluminium	812-335-L	812-335-R
CDS® Hinge 360° plus 24 mm aluminium	812-340-L	812-340-R
CDS® Adapter Hinge 360° plus 16 mm stainless steel	812-350-L	812-350-R
CDS® Adapter Hinge 360° plus 20 mm stainless steel	812-355-L	812-355-R
CDS® Adapter Hinge 360° plus 24 mm stainless steel	812-360-L	812-360-R
CDS® Adapter Hinge 360° plus 20 mm aluminium	812-385-L	812-385-R
CDS® Adapter Hinge 360° plus 24 mm aluminium	812-390-L	812-390-R

Hinges called PLUS include a fine adjustment option. Fine adjustment allows for the infinitely variable adjustment of the redression range. The redression range may also be extended to -15° to enable hyperextension.

The CDS® Nano Hinge

Features

- Individual adjustment of the spring tension up to 2,2 Nm
- Spring tension can be switched on and off without tools and without varying the set spring tension
- The redression range can be adjusted in 10° increments

CDS® Nano Hinges	Art.No. left	Art.No. right
CDS® Nano Hinge 16 mm stainless steel	812-100-L	812-100-R
CDS® Nano Hinge 16 mm aluminium	812-130-L	812-130-R
CDS® Nano Adapter Hinge 16 mm stainless steel	812-150-L	812-150-R
CDS® Nano Adapter Hinge 16 mm aluminium	812-180-L	812-180-R

The ROM Nano Adapter Hinge

Features

- Limits may be adjusted without tools.
- Extension limit can be adjusted between -10° and 100° in 10° steps
- Flexion limits can be adjusted between 0° and 130° in 10° steps
- Immobilisation possible between -10° and 100° in 10° steps

ROM Nano Adapter Hinge	Art.No. left	Art.No. right
ROM Nano Adapter Hinge stainless steel	811-110-S	811-111-S
ROM Nano Adapter Hinge aluminium	811-100-S	811-101-S

The CDS® Micro Hinge

Features

- Individual adjustment of the spring tension up to 0,43 Nm

CDS® Micro Hinge	Art.No.
CDS® Micro Hinge	840-500-S

Materials an Width

	CDS® (3 mm thickness)		CDS® Nano (2 mm thickness)
Aluminium	20 mm	24 mm	16 mm
Stainless steel	16 mm	20 mm	16 mm


Casting options

Hinge

The rods are fixed to the hinges and can be riveted to the cast.

Adapter Hinge

The rods can be unscrewed and removed from the adapter-hinges. So they can be poured into the cast.

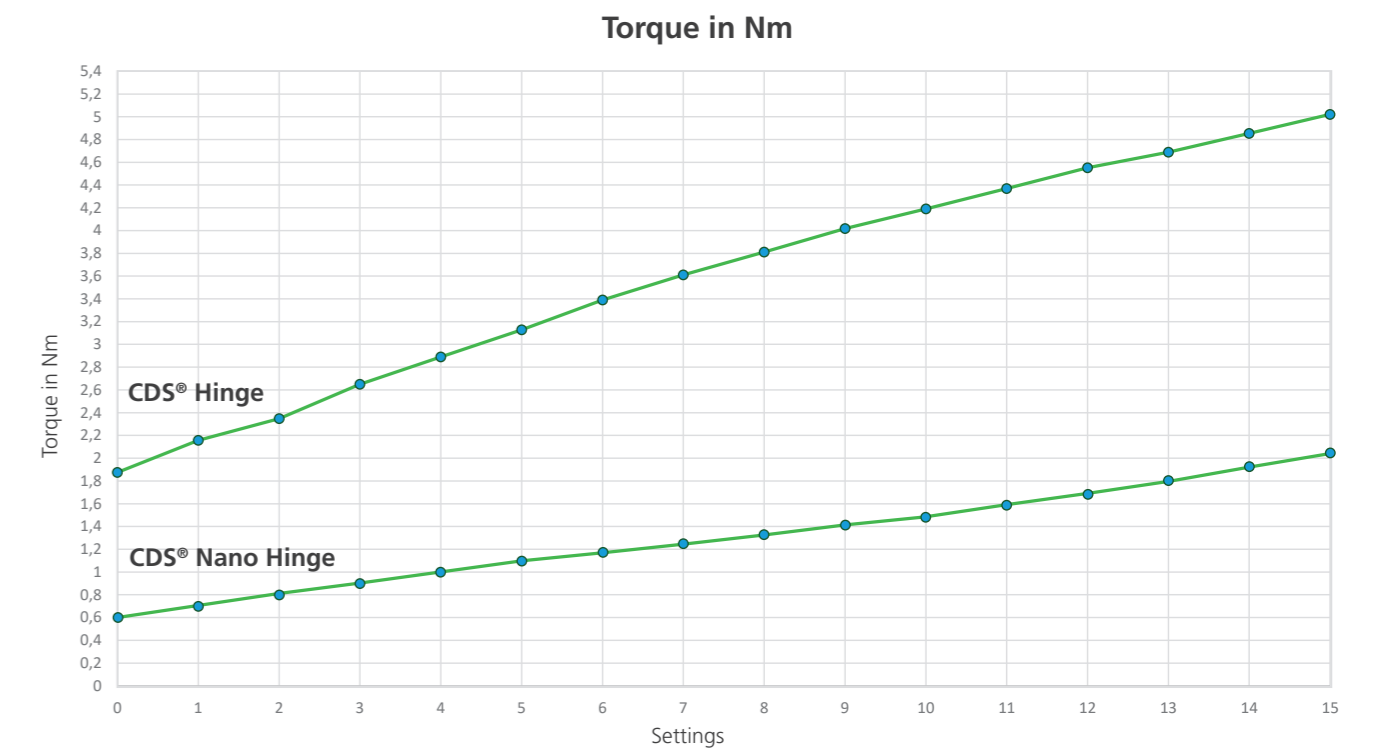

Quality for creativity.



Made in Germany
Hergestellt in Deutschland

en 0220



albrecht GmbH
CHIEMSEESTRASSE 81

D-83233 BERNAU AM CHIEMSEE
GERMANY

Phone + 49 8051 96129 - 0
FAX + 49 8051 96129 - 30

INFO@ALBRECHTGMBH.COM
WWW.ALBRECHTGMBH.COM